


Voltage Sensor

Setup Instructions

Parts Checklist

The Voltage Sensor is only to be used in conjunction with the Data Collector. This is a 0-5V Voltage Probe. To get a valid reading between 0 and 5 volts the red lead must be connected to the side of the circuit with the higher potential. If it is connected in reverse the reading will be near zero.

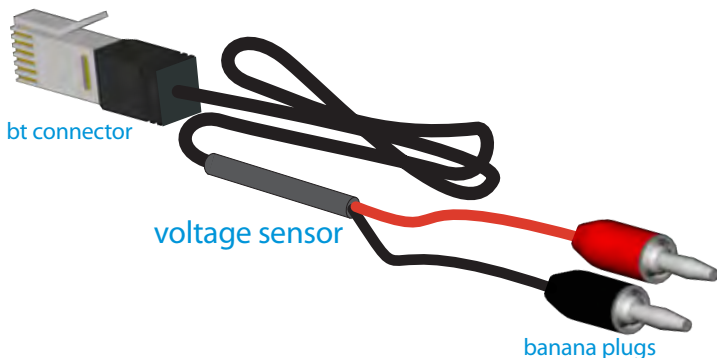

Caution

This experiment deals with electricity. Please do not use the plugs to measure anything over 5 volts. This includes 9V batteries and wall plugs.

Other Activities

For activities, refer to the Physics, A First Course Investigations manual or Teacher's Guide.

Visit us at cposcience.com. If you need technical assistance, please call 1-800-932-5227